

Setrox S Lead

Steroid-Euting, Active Fixation Pacing Endocardial Leads
IS-1 Connector
Technical Manual



CAUTION

Federal (U.S.A.) law restricts this device to sale by, or on the order of,
a physician (or properly licensed practitioner).

© 2015 BIOTRONIK, Inc., all rights reserved.

Contents

1. Device Description	1
2. Indications for Use	3
3. Contraindications	5
4. Warnings and Precautions	7
5. Potential Adverse Events.....	11
6. Clinical Studies	13
7. General Information on Product Handling	15
7.1 Sterilization and Storage	15
7.2 Opening the Sterile Container	15
7.3 Accessories	16
8. Implantation	17
8.1 General Guidelines.....	17
8.2 Measurement of the Capture Threshold and Intracardiac Signals	20
8.3 Lead Connection to the Pulse Generator	21
9. Disclaimer.....	23
10. Technical Specifications	25

1. Device Description

The BIOTRONIK Setrox S transvenous, steroid-eluting, active fixation endocardial lead family is designed for permanent pacing and sensing. Active fixation pacing leads with a bipolar (BP) IS-1 connector configuration are designed for use in conjunction with implantable pulse generators with IS-1 compatible headers. The leads may be used with single- or dual-chamber pulse generators, dual chamber ICDs, CRT-P and D-devices.

Setrox S leads feature an electrically active extendable/retractable fixation helix for use in lead placement. The helix is extended and retracted by rotating the connector pin with a fixation tool. Both the fixation helix and ring electrode are comprised of a platinum/iridium alloy base with a fractal iridium surface. The fractal surface of the lead electrodes creates a larger effective surface area, as a result maximizes the myocardial interface, which is a major factor in determining a lead's sensing characteristics. All leads are multifilar and insulated with medical grade silicone.

The distal tip of the Setrox S lead consists of a steroid-eluting collar, containing 0.75 mg of dexamethasone acetate (DXA). Upon exposure to body fluids, the steroid elutes from the collar into the body tissue by diffusion. Release of the steroid is intended to decrease the inflammatory response at the contact site between the lead tip and the endocardium, thereby decreasing the elevated pacing thresholds of the endocardial lead that often occur after lead implantation.

Setrox S leads have straight distal ends (Setrox S xx) and are intended for placement in either the right atrium or right ventricle. The "xx" represents the lead length in centimeters. The Setrox S leads are available in the following configurations: Setrox S 45, Setrox S 53, and Setrox S 60.

CAUTION

Because of the numerous available 3.2 mm configurations, e.g., the IS-1 and VS-1 standards, lead/pulse generator compatibility should be confirmed with the pulse generator and/or lead manufacturer prior to the implantation of a pacing system.

NOTE:

IS-1, wherever stated in this manual, refers to the international standard, whereby leads and generators from different manufacturers are assured a basic fit. [Reference ISO 5841-3:1992(E)].

See Section 10 for the technical specifications of the Setrox S leads.

2. Indications for Use

BIOTRONIK's Setrox S xx transvenous, steroid-eluting, active fixation endocardial leads are indicated for permanent pacing and sensing. Active fixation pacing leads with a bipolar (BP) IS-1 connector configuration are designed for use in conjunction with implantable pulse generators with IS-1 headers. The leads may be used with single or dual chamber pacing systems, dual chamber ICDs, CRT-P and CRT-D.

The Setrox S xx lead models are intended for placement in either the right atrium or right ventricle.

3. Contraindications

Transvenous endocardial pacing leads are contraindicated in the presence of severe tricuspid valvular disease and in patients with mechanical heart valves. The Setrox S lead is additionally contraindicated for patients who cannot tolerate a single systemic dose of up to 1.3 mg of dexamethasone acetate (DXA).

4. Warnings and Precautions

The performance of a cardiac pacing system depends on proper interaction of its three components: the pulse generator, the lead(s), and the patient. Abnormalities or changes in the electrical properties of any of the three components, or their interfaces with each other, may directly affect function of the entire system. Correct lead implantation is critical to safe and effective performance of the pacing system.

The pacing system may cease to function at any time due to medical and/or technical complications:

Medical Complications

Medical complications of the pacemaker treatment may include, but are not limited to: fibrotic tissue formation, thrombosis, embolism, elevated thresholds, body rejection phenomena, cardiac tamponade, muscle and nerve stimulation, myocardial perforation, erosion of the pulse generator/lead through the skin, infection and pacemaker-induced dysrhythmia (some of which could be life-threatening such as ventricular fibrillation).

Technical Complications

Incorrect operation of the pacing system may be caused by but is not limited to: improper lead placement, lead dislodgement, lead fracture, loss of insulation integrity, battery depletion, or electrical component failure.

Potentially Harmful Therapeutic and Diagnostic Procedures

As an implanted pacing lead is a direct, low resistance path to the myocardium for electrical current, the observance of high standards of electrical safety is required. Electrosurgical instruments, for example, could generate voltages of such amplitude that a direct coupling between the tip of the electrocautery device and the implanted lead may result, possibly inducing myocardial lesions or serious cardiac arrhythmias (e.g., fibrillation).

Some therapeutic and diagnostic procedures (e.g., diathermy, MRI, electrocautery) may result in latent damage to the pacing system. This damage may not be detected when testing the pulse generator function immediately after the procedure, but may become evident at a later time, resulting in pacing system malfunction or failure.

Lithotripsy – Lithotripsy treatment should be avoided since electrical and/or mechanical interference with the pacemaker or ICD is possible. If this procedure must be used, the greatest possible distance from the point of electrical and mechanical strain should be chosen in order to minimize a potential interference with the implant.

Medical Procedures – For any medical procedures that may affect the device (e.g., therapeutic ultrasound, external defibrillation, electrophysiological ablation, HF surgery, lithotripsy), perform a complete follow-up after the procedure.

Therapeutic Ultrasound – Therapeutic ultrasound is not recommended due to possible heating effects of the device at the implant site. If therapeutic ultrasound must be considered, it should not be applied in the immediate vicinity of the implant.

Transcutaneous Electrical Nerve Stimulation (TENS) – Transcutaneous Electrical Nerve Stimulation should be avoided, as it may lead to unintended heart stimulation.

Prevention of Leakage Current Conduction

Pulse generators and testing equipment connected to the lead must be battery-powered. Proper grounding of line-powered devices in the vicinity of the patient is essential to prevent leakage currents arising from such devices to be conducted via the lead's terminal or any other non-insulated part.

Excessive Pressure

Excessive pressure and hyperbaric oxygen therapy should be avoided, as it may cause damage to the implant.

Previously Implanted Leads

It is generally recommended that a chronically implanted endocardial lead not be explanted. If it becomes necessary to abandon a lead, the connector pin should be capped to prevent the transmission of electrical signals to the heart.

Storage Temperature

Storage at temperatures up to 25° C (77° F); excursions permitted from 5° to 55° C (41° to 131° F). Exposure to temperatures outside this range may result in lead malfunction.

Necessary Equipment for Implantation

During implantation the ECG should be recorded; a pacing system analyzer (PSA) and defibrillation equipment should always be readily available.

Handling the Lead

The lead should be handled very carefully at all times. Any severe application of force (bending, stretching, crimping, etc.) may permanently damage the lead. The metal portion of the lead connector should not be touched.

Intrusion of Blood

Avoid intrusion of blood into the lead lumen.

Lead Positioning

If the ICD or pacemaker is implanted underneath the pectoral muscle, ensure that no parts of the lead lie between the ribs and clavicle or between the housing of the implant and the ribs/clavicle. Chafing and pressure on the lead between the housing of the implant and the ribs/clavicle could damage the lead's insulation and thus cause premature failure.

Stylet Compatibility

To ensure compatibility, use only the stylets that are packaged with the lead or approved for use with the lead. Unapproved stylet types may result damage to the lead and/or patient injury.

Stylet Insertion

To avoid damage to the lead, do not insert the stylet too rapidly nor use excessive force when inserting the stylet into the lead.

Lead/Pulse Generator Compatibility

Because of the numerous available 3.2 mm configurations, e.g., the IS-1 and VS-1 standards, lead/pulse generator compatibility should be confirmed with the pulse generator and/or lead manufacturer prior to the implantation of a pacing system.

NOTE:

IS-1, wherever stated in this manual, refers to the international standard, whereby leads and generators from different manufacturers are assured a basic fit. [Reference ISO 5841-3:1992(E)].

Extending/Retracting the Fixation Helix

In the event of previous handling or repositioning of the lead, more than the minimum number of rotations may be necessary to fully extend or retract the helix. Full helix extension should always be verified through fluoroscopy.

If the screw mechanism becomes difficult to handle or if it or the stylet sticks due to repeated extension and retraction of the fixation screw (from repositioning of the lead tip), the lead should be removed and replaced with a new one by following these measures:

- No longer use the screw mechanism.
- Rotate the entire lead with inserted stylet counterclockwise in order to unscrew the lead from the myocardium without using the screw mechanism.

Anchoring Sleeve

Always use an anchoring sleeve (lead fixation sleeve) when implanting a lead. Use of the anchoring sleeve, which is provided with the lead, will lessen the possibility of lead dislodgement and protect the lead body from damage by a securing ligature.

Measuring Intracardiac Signals

Depending on the PSA used, pacing may be interrupted during the measurement of the intracardiac signals.

Chronic Repositioning

It is generally recommended that a chronically implanted endocardial lead not be explanted. Chronic repositioning or removal of active fixation leads may be difficult due to the presence of blood or fibrotic tissue in the helix. If it becomes necessary to remove the lead without successfully retracting the fixation helix, the lead should be rotated counter-clockwise during withdrawal in order to minimize the risk of endothelial laceration. If it becomes necessary to abandon a lead, the connector pin should be capped to prevent the transmission of electrical signals to the heart.

Setscrew Adjustment

The pulse generator's setscrew(s) must be retracted prior to inserting the lead connector. Failure to back off the pulse generator's setscrew(s) may result in damage to the lead(s), and/or difficulty connecting the lead(s).

Cross-Threading Setscrew

To prevent cross-threading the setscrew, do not back the setscrew completely out of the threaded hole. Leave the torque wrench in the slot of the setscrew while the lead is inserted.

Tightening Setscrew

Do not over-tighten the setscrew(s). Use only a torque wrench, which automatically prevents over-tightening.

Sealing Caps

For pacemakers requiring sealing caps, secure a sealing cap over the setscrew(s) to prevent pacemaker malfunction.

5. Potential Adverse Events

Potential complications resulting from the use of endocardial leads include, but are not limited to: thrombosis, embolism, body rejection phenomena, cardiac tamponade, pneumothorax, muscle/nerve stimulation, valve damage, fibrillation, infection, skin erosion and ventricular ectopy. Lead perforation through the myocardium has been rarely observed. The table below summarizes some of the potential symptoms indicating a complication and possible corrective actions:

Table 1: Potential Complications and Corrective Actions

Symptom	Potential Complication	Potential Corrective Action
Loss of pacing or sensing	Lead dislodgement	Reposition lead
	Lead fracture	Replace lead
	Lead insulation defect	Replace lead
	Improper lead / pulse generator connection	Reconnect lead to pulse generator
Increase/ decrease in threshold	Fibrotic tissue formation	Adjust pulse generator output; Replace/reposition lead

6. Clinical Studies

Because the Setrox S lead is based on BIOTRONIK's legally marketed Selox SR lead, engineering and animal tests were performed in lieu of human clinical data.

7. General Information on Product Handling

The following information generally applies to all transvenous implantable leads, but does not attempt to describe all procedures to be followed, precautions to be taken, or contraindications to be considered.

7.1 Sterilization and Storage

The lead and its accessories have been sealed in a double blister plastic container and gas sterilized with ethylene oxide. To assure sterility, the container should be checked for integrity prior to opening. Should a breach of sterility be suspected, return the lead to BIOTRONIK.

The blister package is shipped in an outer box, equipped with a quality control seal and product information label. The label contains the model specifications, technical data, serial number, shelf-life expiration date, sterilization and storage information of the lead.

CAUTION

Storage Temperature

Storage at temperatures up to 25° C (77° F); excursions permitted from 5° to 55° C (41° to 131° F). Exposure to temperatures outside this range may result in lead malfunction.


Handling the Lead

The lead should be handled very carefully at all times. Any severe application of force (bending, stretching, crimping, etc.) may permanently damage the lead. The metal portion of the lead connector should not be touched.


Should a replacement lead be required, contact your local BIOTRONIK representative.

7.2 Opening the Sterile Container

The lead is packaged in two plastic containers, one within the other. The two packages are individually sealed and then sterilized with ethylene oxide gas. Due to the double packaging, the outside of the inner container is sterile and can be removed using standard aseptic technique and placed in the sterile field.

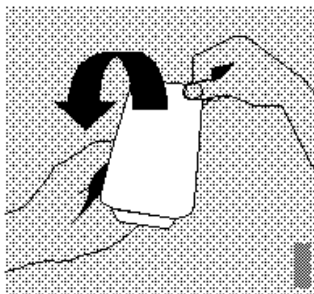

Peel off the sealing paper of the outer unsterile container as indicated by the arrow.

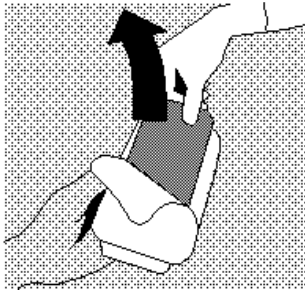

Take out the inner sterile container by the gripping tab and open it by peeling the sealing paper as indicated by the arrow.

7.3 Accessories

Compatible accessories include suture sleeves, stylets, stylet guides, fixation tools and vein lifters.

8. Implantation

The lead tip consists of a silicone rubber collar that contains a steroid agent. The steroid is used to reduce inflammation after the implantation and decrease the postoperative threshold increase caused by inflammation.


8.1 General Guidelines

Setrox S leads come with a regular straight ball-tipped stylet already inserted.

Verifying the Mechanical Function of the Fixation Helix

Prior to implantation, test the operation of the fixation helix:


1. Press both legs of the fixation tool together and place the most distal hole of the fixation tool on the connector pin.

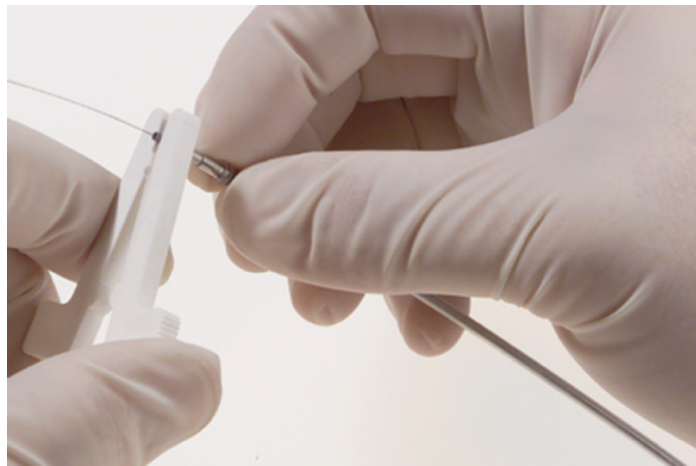

NOTE:

Only employ the screw mechanism using the provided fixation tool, which is clamped to the connector pin. Do not use any other tools or accessories.

2. Practice extending the helix by rotating the tool clockwise until the helix is completely exposed (a minimum of 6-9 complete clockwise rotations of the fixation tool). This number of rotations of the fixation tool should fully extend the helix to approximately 1.8 mm. (See following figure).

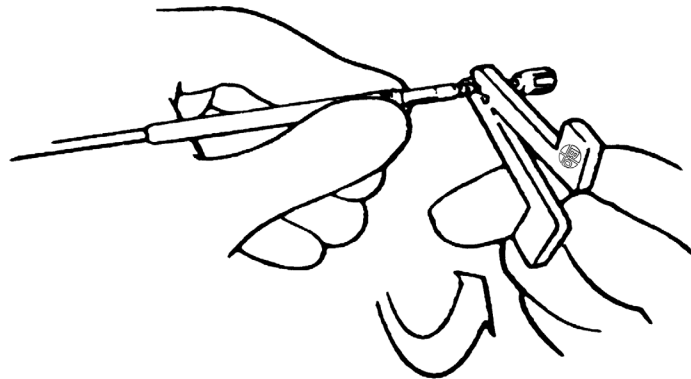

3. Remove the fixation tool from the connector pin and release the proximal end of the lead body. Allow the residual torque in the lead to be relieved.

4. Reattach the fixation tool. Practice retracting the helix by rotating the tool counterclockwise until the helix is retracted in the sheath (a minimum of 6 complete counterclockwise rotations of the fixation tool).

NOTE:

More than the recommended minimum of 6-9 rotations may be necessary to extend or retract the fixation helix.

Venous Access

A number of venous routes are available for implanting an endocardial lead, including the cephalic, subclavian, and external or internal jugular veins. Venous access can be gained by using either the cutdown or subclavian puncture technique.

CAUTION


Lead/Pulse Generator Compatibility

Because of the numerous available 3.2 mm configurations, e.g., the IS-1 and VS-1 standards, lead/pulse generator compatibility should be confirmed with the pulse generator and/or lead manufacturer prior to the implantation of a pacing system.

NOTE:

IS-1, wherever stated in this manual, refers to the international standard, whereby leads and generators from different manufacturers are assured a basic fit. [Reference ISO 5841-3:1992(E)].

5. Position the lead fixation sleeve close to the lead's connector pin.
6. To employ the cephalic cutdown technique, carefully insert the pointed end of the disposable vein lifter into the lumen after opening the vein, and then lift the vein to permit easy lead insertion. (see following figure)

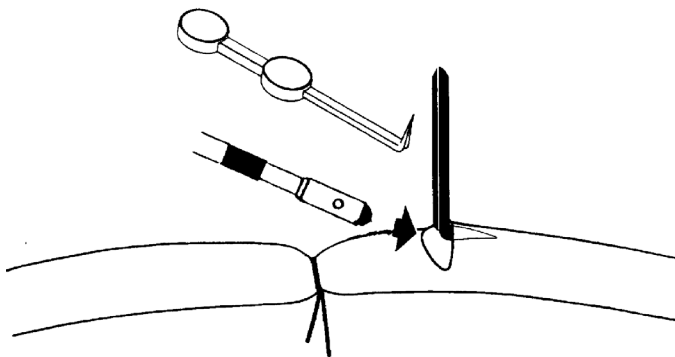

7. If the subclavian route is selected and access by subclavian puncture is preferred, a percutaneous lead introducer should be used. Perform a subclavian puncture and insert a lead introducer into the vein. Remove the introducer guidewire and dilator, leaving the sheath in place in the vein to receive the lead.
8. Advance the stylet all the way into the lead. This will stiffen and straighten the lead, and prevent lead bending.
9. Excessive lead slack can potentially lead to perforation.

NOTE:

Avoid blood contact to the stylet which may make its reinsertion and/or removal difficult. A spare stylet is provided in the inner tray.

NOTE:

BIOTRONIK stylets are designed to be 1 to 2 cm longer than the stylet lumen of the lead body.

CAUTION

Stylet Compatibility

To ensure compatibility, use only the stylets that are packaged with the lead or approved for use with the lead. Unapproved stylet types may result damage to the lead and/or patient injury.

Stylet Insertion

To avoid damage to the lead, do not insert the stylet too rapidly nor use excessive force when inserting the stylet into the lead.

Securing the Electrode into the Endocardium

10. After the electrode has been advanced into a stable position within the atrium or ventricle, leave the stylet in the lead and attach the fixation tool to the connector pin.
11. Achieve permanent lead fixation by rotating the fixation tool clockwise until the helix is completely exposed.

NOTE:

A minimum of 6-9 rotations of the fixation tool are recommended to fully extend (or retract) the fixation helix.

CAUTION


In the event of previous handling or repositioning of the lead, more than the minimum number of rotations may be necessary to fully extend or retract the helix. Full helix extension should always be verified through fluoroscopy.

12. Use fluoroscopy to verify the position of the fixation helix. Both the helix and the tip are visible under fluoroscopy as noted in the figure below.

Helix retracted:



Helix extended:

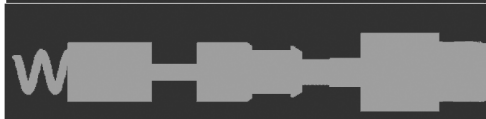


13. Carefully remove the stylet after the electrode has been positioned and fixed into a stable position.

To lessen the possibility of dislodgement, it is recommended to anchor the lead at the incision site where it enters the vein. To facilitate anchoring without damaging the lead insulation or the conductor coil, BIOTRONIK leads are supplied with a silicone rubber fixation sleeve designed with ligature grooves and anchoring tabs.

CAUTION


Anchoring Sleeve

Always use an anchoring sleeve (lead fixation sleeve) when implanting a lead. Use of the anchoring sleeve, which is provided with the lead, will lessen the possibility of lead dislodgement and protect the lead body from damage by a securing ligature.

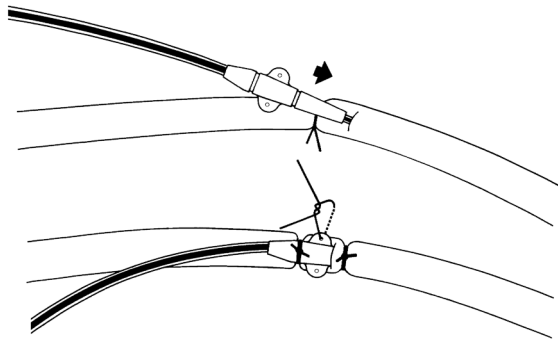

Care should be taken to avoid perforating or otherwise damaging the lead with the stylet, clamps or other surgical instruments. The white marker clamp available separately may be used for identification during a dual lead implantation procedure.

The notch at the side of the stylet guide permits connection of a surgical cable alligator clip to the lead connector pin for measurement of the capture threshold and the intracardiac signal. After this, the pacemaker is connected to the lead as described on page 28 and in the appropriate technical manual for pacemaker implantation.

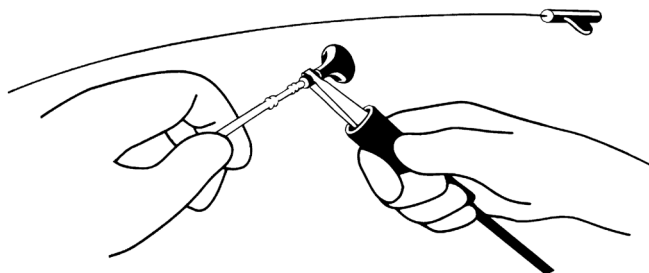

NOTE:

Special care should be taken not to damage the sealing rings of the lead connector, particularly when making contact with the ring electrode of a bipolar lead with an IS-1 connector (letter designation "BP") by means of an alligator clip.

8.2 Measurement of the Capture Threshold and Intracardiac Signals

Low capture thresholds and adequate sensing of intracardiac signal amplitudes indicate appropriate lead placement. The BIOTRONIK Pacing System Analyzer or a comparable device may be used to accurately determine pacing capture thresholds, as well as intracardiac signals, in accordance with implant characteristics of implantable pulse generators.

WARNING

Prevention of Leakage Current Conduction

Pulse generators and testing equipment connected to the lead must be battery-powered. Proper grounding of line-powered devices in the vicinity of the patient is essential to prevent leakage currents arising from such devices to be conducted via the lead's terminal or any other non-insulated part.

Previously Implanted Leads

It is generally recommended that a chronically implanted endocardial lead not be explanted. If it becomes necessary to abandon a lead, the connector pin should be capped to prevent the transmission of electrical signals to the heart.

Generally, the electrode position is regarded acceptable if the values for capture and sensing threshold meet the following industry-recommended values:

Measurement	Atrial	Ventricular
Maximum acute capture threshold at pulse duration setting ≤ 0.5 ms	1.5 V	1.0 V
Minimum acute sensing amplitudes	2.0 mV	5.0 mV
Impedance range (0.5 ms, 4.8 V)	500 to 2000 Ω	

8.3 Lead Connection to the Pulse Generator

1. Withdraw the stylet and stylet guide before connecting the lead to the pulse generator.
2. Using the pulse generator torque wrench, retract the setscrew(s) counterclockwise far enough to ensure unimpeded insertion of the lead connector(s) into the port(s).

CAUTION

Setscrew Adjustment

The pulse generator's setscrew(s) must be retracted prior to inserting the lead connector. Failure to back off the pulse generator's setscrew(s) may result in damage to the lead(s), and/or difficulty connecting the lead(s).

3. Insert the lead connector into the pulse generator's header receptacle following the manufacturer's directions for lead insertion. The method of insertion depends upon the header configuration of the pulse generator.
4. Securely tighten the setscrew(s) of the pulse generator connector.

CAUTION

Tightening Setscrew

Do not over-tighten the setscrew(s). Use only a torque wrench, which automatically prevents over-tightening.

5. Assure correct insertion of the connector pin by verifying that the tip of the pin is visible past the header connector block.
6. Visually and mechanically check the integrity of the connection after tightening the setscrew(s).



CAUTION

Sealing Caps


For pacemakers requiring sealing caps, secure a sealing cap over the setscrew(s) to prevent pacemaker malfunction.

9. Disclaimer

BIOTRONIK leads, lead extensions, adapters and accessories used in connection with these devices (referred to as: leads and accessories) have been qualified, manufactured and tested in accordance with well-proven and accepted standards and procedures. The physician should be aware, however, that leads and accessories may be easily damaged by improper handling or use. Except as set forth in BIOTRONIK's Lead Limited Warranty, BIOTRONIK makes no express or implied warranties for its leads and accessories.

10. Technical Specifications

Setrox S Steroid-Eluting, Bipolar Implantable Endocardial Leads with Active Fixation

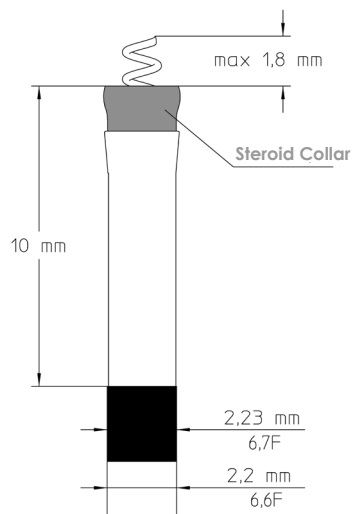

Setrox S 45
Setrox S 53
Setrox S 60

Model (xx = length in cm)	Order Number
Setrox S	
Setrox S 45	350 973
Setrox S 53	350 974
Setrox S 60	350 975

Setrox S xx Technical Parameters	
Connection	IS-1
Polarity	Bipolar
Configuration	Straight
Fixation Helix	
Fixation	Electrically active helix, retractable
Extended	max. 1.8 mm
Surface Area	4.5 mm ²
Material	70% Pt / 30% Ir
Surface structure	Iridium, Fractal
Steroid	
Steroid	Dexamethasone acetate (DXA)
Steroid Amount	0.75 mg
Steroid Carrier	Silicone Rubber
Ring Electrode	
Surface area	17.5 mm ²
Material	90% Pt / 10% Ir
Surface structure	Iridium, Fractal
Tip to Ring Distance	10 mm
Lead Body	
Insulation	Silicone Rubber
Conductor ¹	MP35N
Lead Body diameter	2.2 mm (6.6 F)
Maximum diameter	2.23 mm (6.7 F)
Number of wire filaments	4
Lengths	45 / 53 / 60 cm
Wire Resistance (Tip) Lead Length = 45 cm	0.65 Ω/cm
Wire Resistance (Ring)	2.04 Ω/cm
Recommended Lead Introducer	7 F ²

(Footnotes)

1 MP35N[®] is a trademark for special cobalt-chrome-nickel alloys.

2 Use of a 7 French introducer is recommended if a guide wire is not used during the implantation procedure.

Setrox S

Technical Manual



Manufactured by:

BIOTRONIK SE & Co. KG
Woermannkehre 1
12359 Berlin Germany

M4113-E 09/15
© 2015 BIOTRONIK, Inc. All rights reserved.
MN039r1 9/17/15

BIOTRONIK, Inc.
6024 Jean Road
Lake Oswego, OR 97035-5369
(800) 547-0394 (24-hour)
(800) 291-0470 (fax)
www.biotronik.com

 **BIOTRONIK**
excellence for life